



Calentador de sueros EC230*l* Enthermics (U.S.A.)

- ✓ Calentador de sueros
- ✓ Control de Tº electrónico
- ✓ Capacidad 16 Lt capacidad
- ✓ WarmRight® entrega la temperatura correcta para el producto correcto para fluidos de irrigación o inyección
- ✓ Alarmas visuales y audibles
- ✓ Control de T^o electrónico
- ✓ Sistema de convección forzada
- ✓ Rango Tº fluidos INJ 37-40°C IRR 37-66°C
- ✓ Sensibilidad +0/1,1°C
- ✓ Dimensiones
- √ 409,0 mm ancho
- ✓ 707,9 mm fondo
- √ 492,8 mm alto

Peso 38 kg.







$\textbf{EC230} \ell$ Fluid Warming Cabinet

16 liter fluid chamber capacity
Glass door

WarmRight® provides the right temperature for the right product: injection or irrigation fluids

Injection (INJ) or irrigation (IRR) fluid warming modes 37-40°C (98-104°F) adjustable INJ mode 37-66°C (98-150°F) adjustable IRR mode



Specifications:

The single-chambered warming cabinet is constructed with 20 gauge stainless steel exterior casing and door. Window door allows observation of inventory with the door closed. The cabinet is warmed using a low heat-density electrothermal cable array. The electronic control ensures an accurate chamber temperature within +0/-1.1°C (+0/-2°F) of the set point for temperatures set between 37 - 43°C $(98 - 110^{\circ}F)$ and within $+0/-1.7^{\circ}C$ $(+0/-3^{\circ}F)$ of the set point for temperatures set between 43 - 66°C (110 - 150°F). The cabinet is equipped with one (1) white, epoxy-coated wire basket with a 16 liter capacity to accommodate most commonly used surgical fluids packaged in bags or bottles. The cabinet is furnished with four (4) 1-1/4" (31mm) non-skid rubber feet.

Special Information:

- Right-hand door swing is standard. Lefthand door is available as a special order.
- Stacking combinations listed under Options and Accessories. When ordering a combination of units to stack order one set of stacking hardware.
- WarmWatch™ available as special order.

Control:

The chamber is controlled by one (1) electronic control which consists of a four digit L.E.D. display, on/off key, heat indicator L.E.D, increase and decrease keys, built-in lock-out feature and a series of prompt sequence indicators. Individual temperature controls for irrigation fluids (IRR) and injection fluids (INJ) allow for precise temperature settings and can be easily set to operate in Celsius or Fahrenheit. The IRR temperature may be adjusted from 37 - 66°C (98 - 150°F) and the INJ temperature can be adjusted from 37 - 40°C (98 - 104°F). An alarm will sound if the temperature exceeds 6°C (10°F) over the set-point temperature. In the event of a power failure, the cabinet will remember its programming and begin to operate as before when power is restored. A warming shut-off system, separate from the electronic control, is an additional safety feature.

WarmWatch™ option:

This optional package includes a *dual temperature readout* that provides visual verification of the unit's set-point and actual temperatures; *reporting capabilities* to record and report warmer performance data; and a *USB port* on the control panel for easy data download.

Interior Capacity


EC230ℓ Fluid Warming Cabinet

WarmWatch features:

- An independent temperature sensing system with separate RTD probe, in addition to the standard RTD control probe.
- A dual temperature display that shows set-point temperature and the independently verified, actual chamber logged temperature.
- An electronic control with dedicated microprocessor for PID temperature control and verification temperature sensing systems; continuous electronic monitoring of machine functions and accuracy analysis; and data logging.
- A flush-mounted USB port on the control panel that writes data from memory to a USB flash drive for transfer to a computer. (Computer not included.) Connect USB flash drive to a computer to retrieve written reports and summary.

Contact your Enthermics sales representative for more information.

16 liters

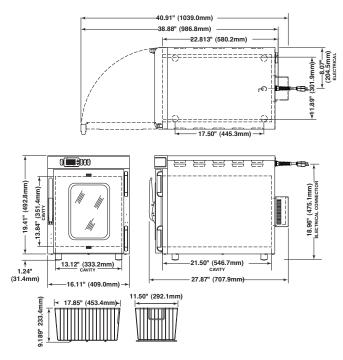

Net Weight	84 lb (38 kg)			
Ship Weight	135	g)		
I	Electrical			
120 V.A.C. 50/60 Hz, 1ph 650 Watts, 5.4 Amps	1 1	NEMA 5-15P 15A - 125V Plug Hospital Grade		, Will the
230 V.A.C. 50/60 Hz, 1ph 800 Watts, 3.5 Amps		BS 1363 Plug (UK only) CEE 7/7 220-230V Plug		C€
Options	and Accesso	ries		
Cart with Casters - Height 41-1/2" (1054mm) EC230 <i>l</i> and EC230 units fit on top of cart and in lower alcove			Contact factory	
Casters*, 3" (76mm)			14227	
Casters*, 5" (127mm)			4007	
Combination Lock Kit			5008361	
Leg Kit, 6" (152mm)			5205	
Lock for Door Handle			LK-2763	
Stacking Hardware EC230 over EC230l over EC230l			5012145	
*Platform (E5089) must be ordered as	a base when order	ing caster	rs or leas for th	nis uni

Capacity and Weights







www.enthermics.com

An ISO 13485:2003 certified company

W164 N9221 Water Street • PO Box 443 Menomonee Falls, WI 53052-0443 Phone 262.251.8356 • 800.862.9276 Fax 262.251.7067

Email generalinfo@enthermics.com